


TempTale® Specifications and Use



DRY ICE FLEX PROBE
#N00006244
-60°C to -90°C
(Initial Pfizer shipments only)



TempTale®4 Reader
Downloads time-and-temperature data
recorded by the monitor to your PC.

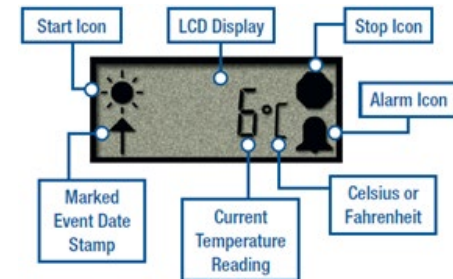

WHITE LABEL

-30°C to 70°C

USAMMA-DOC uses this multi-use model.
(DLA-TSM uses the single-use, colored label models to the right)

ORANGE LABEL

#T4100-01-360

2°C to 8°C

(Refrigerated)
(ex. Flu Vaccines/Novavax/etc.)

YELLOW LABEL

#T4400-02-360

-20°C to 8°C

(Hybrid)
(ex. MMR)

LIGHT BLUE LABEL

#T4400-03-360

15°C to 30°C

(Controlled room temperature)
(ex. Antivirals)

PURPLE LABEL

#T4400-05-360

-25°C to -15°C

(Moderna/Pfizer Redistributions)
Mostly used for OCONUS locations.

- TempTales® validate packaged temperatures throughout the shipping process and provide pertinent data to both customer and shipper to ensure material viability.
- Immediately upon receipt (while still in box), follow DLA-TSM/USAMMA-DOC instructions provided inside vaccine shipping containers for stopping and returning TempTales®.
- Note: the TempTales® that are used for influenza vaccines are the white label and the orange label.**
- Note: for redistributions, ensure the TempTale® you are using is the right temperature range for the specific vaccine(s)/temperature sensitive medical products (TSMP) being shipped.
- Most DoD-approved TempTales® are the green-cased TempTale®4 versions (black cases should no longer be used)
 - Manufacturer issued calibration certificate / 30 day max recording / 1 year battery life / Unit cost: approx. \$100

Manufacturer:

Sensitech Inc.
800 Cummings Center
Beverly, MA 01915

Client Service Associate
Tel: 978-720-2664
Fax: 978-524-4018

(Part #s above are for ordering at Gov't rate.)